

MIPAR™

Network License Activation Instructions


CLIENT SIDE

License Keys may be activated directly within MIPAR using the guided Activation Window.

This window appears after launching MIPAR if there are no active licenses. You may also access it in any app from Help > Activate License.

1 Allow Access to Server

- Go to your License Manager's [Configuration Page](#)
- Click the **Access to Remote License Managers** tab
- Check **Allow Access to Remote Licenses**

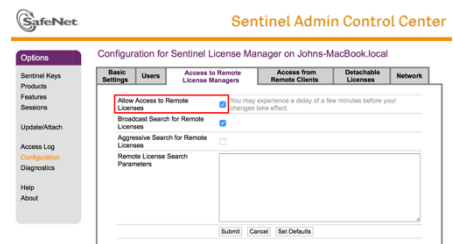

2 Choose Search Type

- Enter **Remote License Search Parameters**
Defines specific machines to search – must be used if access restrictions were set on server

Specify data as IPv4 addresses (for example, 10.1.1.17), Broadcast addresses (for example, 10.1.1.255), or machine names (for example, hklm001.ecomp.com)

Retrieve IPv4 address on server machine using "ipconfig" command in Command Prompt.

Each entry must be on a separate line.

- Click **Submit** to apply any changes

